

"trace & catch"

Instructions for Use

SHIKARI[®] (S-ATD)

Qualitative Antibodies to Denosumab **FIISA**

Enzyme immunoassay for determination of qualitative antibodies to denosumab in serum and plasma samples

REF	DEN-QLS-PRO		
Σ	96 tests		
X	Shipment 10-30°C, Store 2-8°C		
***	MATRIKS BIOTECHNOLOGY CO., LTD. Gazi Universitesi Teknoplaza CZ-26/27 06830 Golbasi Ankara / TURKEY Tel +90 312 485 42 94 info@matriksbiotek.com		
Œ	IVD 🎯 🕈 🔆 🚺		

1. Intended Use

 ${\sf SHIKAR}|^{\otimes}$ Qualitative Antibodies to Denosumab ELISA has been especially developed for the qualitative analysis of antibodies to denosumab in serum and plasma samples. SHIKARI® Qualitative Antibodies to Denosumab ELISA is optimized with Prolia®.

2. General Information

Denosumab is a novel, fully human IoG2 monoclonal antibody specific to receptor activator of nuclear factor kappa-B ligand (RANKL), suppresses bone resorption markers in patients with a variety of metastatic tumours and is being investigated in multiple clinical trials for the and treatment of bone metastases. Chemically, it consists of 2 heavy and 2 light chains. Each light chain consists of 215 amino acids. Each heavy chain consists of 448 amino acids with 4 intramolecular disulphides

Denosumab is designed to target RANKL (RANK ligand), a protein that acts as the primary signal to promote bone removal/resorption. In many bone losses, RANKL overwhelms the body's natural defence against bone destruction. Denosumab prevents RANKL from activating its receptor, RANK, on the surface of osteoclasts and their precursors, Prevention of the RANKL/RANK interaction inhibits osteoclast formation, function, and survival, thereby decreasing bone resorption and increasing bone mass and strength in both cortical and trabecular bone

Therapeutic drug monitoring (TDM) is the clinical practice of measuring specific drugs at designated intervals to maintain a constant concentration in a patient's bloodstream, thereby optimizing individual dosage regimens. The indications for drug monitoring include efficacy, compliance, drug-drug interactions, toxicity avoidance, and therapy cessation monitoring. Additionally, TDM can help to identify problems with medication compliance among noncompliant patient cases.

Biologic medicinal products (biologics) have transformed treatment landscapes worldwide for patients with haematological or solid malignancies with the 21st century. Today, as data exclusivity periods of first wave biologics approach expiration/have expired, several biosimilar products (i.e., biologics that are considered to be similar in terms of quality, safety and efficacy to an approved 'reference' biologic) are being developed or have already been approved for human use

Like all biologics, biosimilars are structurally complex proteins that are typically manufactured using genetically engineered animal, bacterial or plant cell culture systems. As a consequence of this molecular complexity and the proprietary nature of the manufacturing process, which will inevitably result in the use of different host cell lines and expression systems as well as related differences in manufacturing conditions, it is not possible to manufacture exact copies of a reference biologic.

When administered to patients, all therapeutic proteins have the potential to induce an nwanted immune response (i.e., to stimulate the formation of antidrug antibodies [ADAs]). The impact of immune responses can range from no apparent effect to changes in pharmacokinetics, loss of effect and serious adverse events. Furthermore, the immunogenicity profile of a biologic can be significantly altered by even small differences in its manufacturing process that are accompanied by a change in product attributes, as well as differences in dosing schedules, administration routes or patient populations.

SHIKARI® ELISA kits can be used for drug level and anti-drug antibodies measurements. SHIKARI[®] denosumab ELISA products:

Brand	Description		Product Code	
SHIKARI [®] (Q-DEN)	Denosumab	Free drug	DEN-FD-PRO	
SHIKARI [®] (S-ATD)	Denosumab	Antibody screening - Qualitative	DEN-QLS-PRO	
SHIKARI [®] (S-ATD)	Denosumab	Antibody screening - Quantitative	DEN-QNS-PRO	

Check the web page for the whole product list www.matriksbiotek.com

3. Test Principle

Solid phase enzyme-linked immunosorbent assay (ELISA) based on the sandwich principle. Controls and samples (serum or plasma) are incubated in the microtiter plate coated with the drug denosumab. After incubation, the wells are washed. Then, horse radish peroxidase (HRP) conjugated probe is added and binds to denosumab antibodies captured by the drug denosumab on the surface of the wells. Following incubation wells are washed and the bound enzymatic activity is detected by addition of tetramethylbenzidine (TMB) chromogen substrate. Finally, the reaction is terminated with an acidic stop solution. The colour developed is proportional to the amount of denosumab antibodies in the sample or controls. The results can be evaluated with using cut-off value.

4. Warnings and Precautions

For professional use only

- In case of severe damage of the kit package please contact Matriks Biotek® or your supplier in written form, latest one week after receiving the kit. Do not use damaged components in test runs but keep safe for complaint related issues.
- Obey lot number and expiry date. Do not mix reagents of different lots. Do not use expired
- Before starting the assay, read the instructions completely and carefully. Use the valid version of the package insert provided with the kit. Be sure that everything is understood. For further information (clinical background, test performance, automation protocols, alternative applications, literature, etc.) please refer to the local distributor.
- Follow good laboratory practice and safety guidelines. Wear lab coats, disposable latex gloves and protective glasses where necessary.
- All reagents of this kit containing human serum or plasma (standards etc.) have been tested and were found negative for HIV I/II, HBsAg and Anti-HCV. However, a presence of these or other infectious agents cannot be excluded absolutely and therefore reagents should be treated as potential biohazards in use and for disposal
- Reagents of this kit containing hazardous material may cause eve and skin irritations. See "Materials supplied", SDS and labels for details.
- Chemicals and prepared or used reagents must be treated as hazardous waste according the national biohazard safety guidelines or regulations.

5. Storage and Stability

The kit is shipped at ambient temperature (10-30°C) and should be stored at 2-8°C for long term storage. Keep away from heat or direct sunlight. The strips of microtiter plate are stable up to the expiry date of the kit in the broken, but tightly closed bag when stored at 2-8°C.

6. Specimen (Collection and Storage)

Serum, Plasma (EDTA, Heparin)

The usual precautions for venipuncture should be observed. Do not use grossly baemolytic icteric or lipemic specimens. Samples appearing turbid should be centrifuged before testing to remove any particulate material. Avoid repeated freeze-thaw cycles for serum/plasma

Samples should be diluted with the dilution rate given in the "Pre-test setup instructions before the test.

Drug infusions may camouflages/mask the presence of antibody to drugs in serum/plasma samples. Therefore, blood sampling time is critical for detection of antibodies. It is recommended to take the blood sample just before the scheduled dose (trough specimen)

Storage	2-8°C	-20°C
Stability (serum/plasma)	2 days	6 months

7.	Materials Supplied
----	--------------------

Microtiter Plate	1 x 12 x 8	Microtiter plate Break apart strips. Microtiter plate with 12 rows each of 8 wells coated with denosumab.	
Controls	1 mL (negative) 1 mL (positive)	Control negative and positive Ready to use. Contains human serum and stabilizer, <0,1% NaN ₃ .	
Assay Buffer	1 x 50 mL	Assay buffer Ready to use. Blue coloured. Contains proteins, <0,1% NaN ₃ .	
Conjugate	1 x 12 mL	Horse radish peroxidase conjugated probe Ready to use. Red coloured. Contains HRP conjugated probe, stabilizer and preservatives.	
Substrate	1 x 12 mL	TMB substrate solution Ready to use: Contains 3,3',5,5'- Tetramethylbenzidine (TMB).	
Stop Buffer	1 x 12 mL	TMB stop solution Ready to use. 1N HCI.	
Wash Buffer	1 x 50 mL	Wash buffer (20x) Prepared concentrated (20x) and should be diluted with the dilution rate given in the "Pre- test setup instructions" before the test. Contains buffer with tween 20.	
Foil	2 x 1	Adhesive Foil For covering microtiter plate during incubation	

8. Materials Required but Not Supplied

Micropipettes and tips

Calibrated measures

- Tubes for sample dilution

- Wash bottle, automated or semi-automated microtiter plate washing system

- Microtiter plate reader capable of measuring optical density with a photometer at OD 450nm with reference wavelength 650 nm (450/650 nm)

- Distilled or deionised water, paper towels, pipette tips and timer

9. Procedure Notes

- Any improper handling of samples or modification of the test procedure may influence the results. The indicated pipetting volumes, incubation times, temperatures and pretreatment steps must be performed strictly according to the instructions. Use calibrated pipettes and devices only

 Once the test has been started, all steps should be completed without interruption. Make sure that required reagents, materials and devices are prepared ready at the appropriate time. Allow all reagents and specimens to reach room temperature (18-25°C) and gently swirl each vial of liquid reagent and sample before use. Mix reagents without foaming.

 Avoid contamination of reagents, pipettes and wells/tubes. Use new disposable plastic pipette tips for each reagent, standard or specimen. Do not interchange caps. Always cap not used vials. Do not reuse wells/tubes or reagents.

Use a pipetting scheme to verify an appropriate plate layou

Incubation time affects results. All wells should be handled in the same order and time sequences. It is recommended to use an eight-channel micropipette for pipetting of solutions in all wells.

 Microplate washing is important. Improperly washed wells will give erroneous results. It is recommended to use a multichannel pipette or an automatic microplate washing system. Do not allow the wells to dry between incubations. Do not scratch coated wells during rinsing and aspiration. Rinse and fill all reagents with care. While rinsing, check that all wells are filled precisely with wash buffer, and that there are no residues in the

Humidity affects the coated wells/tubes. Do not open the pouch until it reaches room temperature. Unused wells/tubes should be returned immediately to the resealed pouch including the desiccant.

10. Pre-test Setup Instructions

Preparation of components

	-	
Component	Wash buffer (Must be prepared before starting assay procedure)	
Dilute	10 mL (e.g.)	
With	Up to 200 mL	
Diluent	Distilled water	
Dilution Ratio	1/20	
Remarks	Warm up 37°C to dissolve crystals. Mix vigorously.	
Storage	2-8°C	
Stability	2 weeks	


- Dilution of samples

Sample	Serum/Plasma	
Diluent	Assay buffer	
Dilution Ratio	1/10	
Remarks 1/10 dilution 20 µL sample + 180 µL assay buffer		

Patient samples with a concentration of drug above the measuring range are to be rated as > "Highest Standard (Standard A)". The result must not be extrapolated. The patient sample in question should be further diluted with assay buffer and retested.

11. Test Procedure

	Total assay time: 140 minutes		
1	Pipette 100 μL of each "Negative control", "Positive control" and diluted samples into the respective wells of microtiter plate		
	Wells A1: Negative control* B1: Negative control* C1: Positive control D1 and on: Samples *It is advised to run more than one "Negative control" samples. Negative control studies can be duplicated or triplicated in order to take the mean value		
2	Cover the plate with adhesive foil		
	Briefly mix contents by gently shaking the plate		
	Incubate 60 minutes at room temperature (18-25°C)		
3	Remove adhesive foil		
	Discard incubation solution		
	Wash plate three times each with 300 µL "Wash Buffer"		
	Remove excess solution by tapping the inverted plate on a paper towel		
4	Pipette 100 μL "Conjugate" into each well		

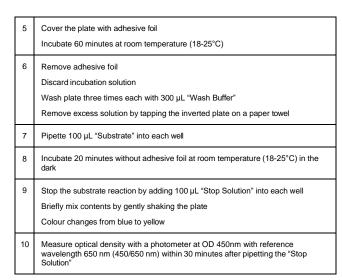

12. Quality Control

The test results are only valid if the test has been performed following the instructions. Moreover, the user must strictly adhere to the rules of GLP (Good Laboratory Practice) or other applicable standards/laws. For the run to be valid, the OD 450/650 nm of positive control should be >1,000 and the OD 450/650 nm of each negative control should be <0,200. In case of any deviation the following technical issues (but not limited to) should be reviewed: Expiration dates of reagents, storage conditions, pipettes, devices, incubation conditions, washing methods, etc.

13. Calculation and Interpretation of Results

The results are evaluated by a cut-off value which is estimated by multiplying the mean OD 450/650 nm of the negative controls by 3.

```
e.g.
```

If "Sample OD 450/650 / the mean negative control OD 450/650 ≥3"

then the sample is $\ensuremath{\mathsf{POSITIVE}}$

If "Sample OD 450/650 / the mean negative control OD 450/650 <3"

then the sample is NEGATIVE

Note: The cut-off information provided with this kit can only be considered as a recommendation. Cut-off values must be calculated/set or verified according to scientific standards by the users/laboratories.

14. Analytical Performance

- Specificity: There is no cross reaction with native serum immunoglobulin
- Precision: Intra-assay and inter-assay CVs <30%
- Cut-off: Cut-off values must be calculated/set or verified according to scientific standards by the users/laboratories.

The "Quality control certificate" contains lot specific analytical performance data and is supplied separately with each kit. If some further analytical performance data is needed, please refer to the local distributor.

15. Automation

 $\mathsf{SHIKARI}^{\textcircled{O}}$ Qualitative Antibodies to Denosumab ELISA is also suitable to run on automated ELISA processors.

16. Symbols and Cautions

	Manufacturer	X	Temperature limitation
~	Production date		See instruction for use
><	Expiry date	A.	Caution
LOT	Lot number	IVD	In vitro diagnostic medical device
REF	Catalog number	CONTROL	Control
۲	Do not use if package is damaged	CONTROL -	Negative control
类	Keep away from sunlight	CONTROL +	Positive control
Ť	Keep dry	∇	Number of tests

Cautions: The performance of the kit can be achieved by fully complying with the instructions. Modifications on the test procedure can affect the results and these kinds of changes will not be charged as regular complaints. This product is for professional use only and must be used for "Intended use" that is given in the instructions for use. The results themselves should not be the only reason for any therapeutically consequences. They must be calculated/set according to scientific standards by the users/laboratories. Information in the instructions about cut-off, etc. performance characteristics, can only be considered as a recommendation and does not give any responsibility to the manufacturer.

Limitations of liability: The manufacturer's liability is limited to the purchase price of the product in all circumstances. The manufacturer cannot be held responsible for damage to the patient, lost profit, lost sales, damage to property or any other incidental or consequential loss.

Technical support and complaints: Technical support can be given upon request. If there is a problem with the product, complaints must be sent written to <u>info@matriksbiotek.com</u> with the technical data (if available) like standard curve, control results, etc. After the necessary examination, written reply will be given.

17. References

Amgen Prolia® (denosumab) product insert.

 Malan J, Ettinger K, Naumann E, Beirne OR. The relationship of denosumab pharmacology and osteonecrosis of the jaws. Oral Surg Oral Med Oral Pathol Oral Radiol 2012 Dec;114(6):671-6.

18. Revision summary

Revision no	Release date	Explanation
01	03.02.2020	New documentation
02	14.10.2021	Change Test Procedure
03	22.12.2022	Sections 7, 10, 11 and 12 have been revised Document code has been changed



 Version
 IFU-0111
 Revision No / Date
 03/22.12.2022